

# Anti-ATM Antibody (R3S31)


## Summary

<b>Catalog No.</b>	RHG48103
<b>Clone ID</b>	R3S31
<b>Host species</b>	Mouse
<b>Tested applications</b>	ELISA: 1:10000, FCM: 1:200-1:400, IF: 1:200-1:1000, IHC: 1:200-1:1000
<b>Species reactivity</b>	Human, Mouse, Rat, Monkey
<b>Form</b>	Liquid
<b>Storage buffer</b>	0.01M PBS, pH 7.4, 0.05% Sodium Azide.
<b>Concentration</b>	1 mg/ml
<b>Purity</b>	>95% as determined by SDS-PAGE.
<b>Clonality</b>	Monoclonal
<b>Isotype</b>	IgG2b
<b>Applications</b>	ELISA, FCM, IF, IHC
<b>Target</b>	A-T mutated, Ataxia telangiectasia mutated, Serine-protein kinase ATM, ATM
<b>Purification</b>	Protein A/G purified from cell culture supernatant.
<b>Endotoxin level</b>	Please contact with the lab for this information.
<b>Accession</b>	Q13315
<b>Stability and Storage</b>	Use a manual defrost freezer and avoid repeated freeze-thaw cycles. Store at 4°C short term (1-2 weeks). Store at -20°C 12 months. Store at -80°C long term.


**Note**

For research use only.


## Data Image

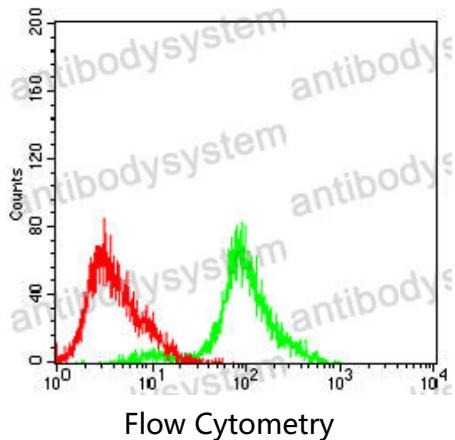

Flow cytometric analysis of Hela cells using ATM mouse mAb (green) and negative control (red).

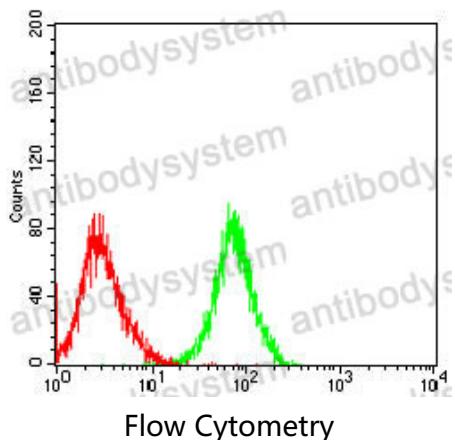


Flow cytometric analysis of COS-7 cells using ATM mouse mAb (green) and negative control (red).

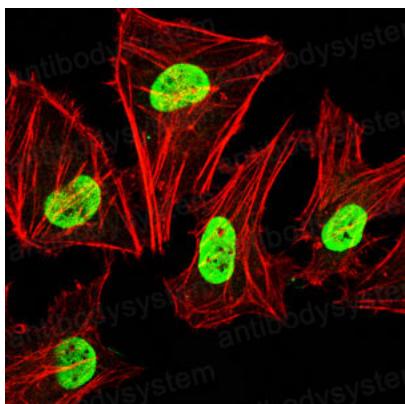


Immunofluorescence

Immunofluorescence analysis of Hela cells using ATM mouse mAb (green). Blue: DRAQ5 fluorescent DNA dye. Red: Actin filaments have been labeled with Alexa Fluor- 555 phalloidin.

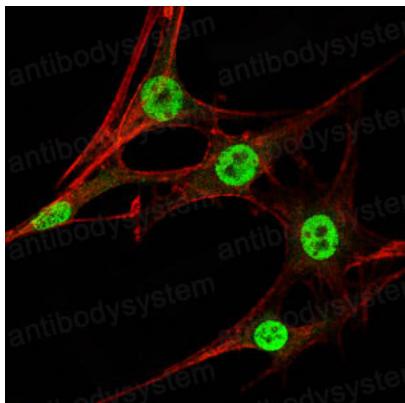


Immunofluorescence

Immunofluorescence analysis of NIH/3T3 cells using ATM mouse mAb (green). Blue: DRAQ5 fluorescent DNA dye. Red: Actin filaments have been labeled with Alexa Fluor- 555 phalloidin.

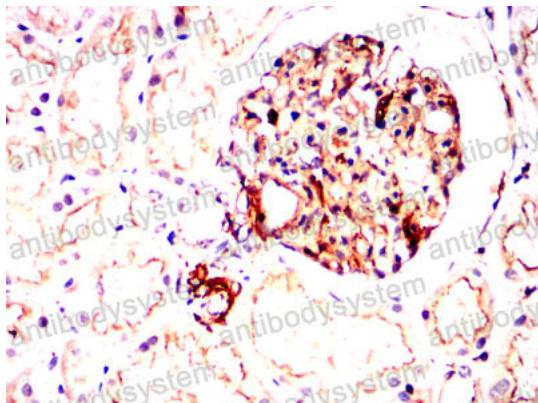


Immunohistochemical

Immunohistochemical analysis of paraffin-embedded human kidney tissues using ATM mouse mAb with DAB staining.

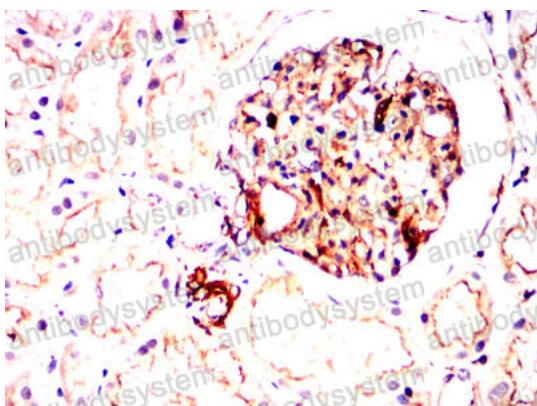


Immunohistochemical

Immunohistochemical analysis of paraffin-embedded human kidney tissues using ATM mouse mAb with DAB staining.

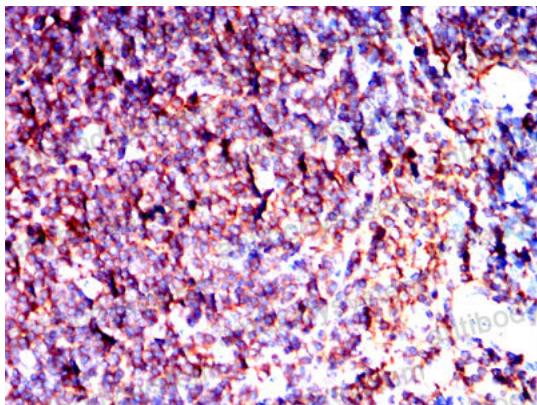


Immunohistochemical

Immunohistochemical analysis of paraffin-embedded human gastric cancer tissues using ATM mouse mAb with DAB staining.

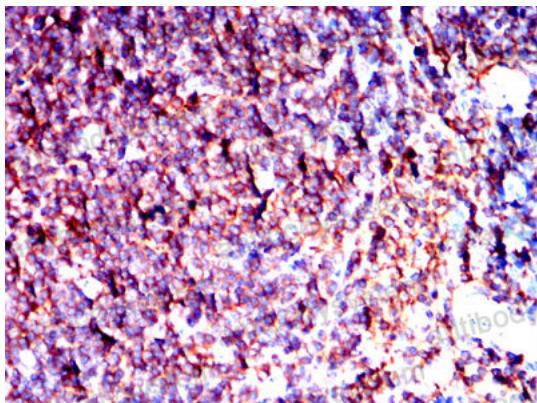


Immunohistochemical

Immunohistochemical analysis of paraffin-embedded Mouse spleen tissues using ATM mouse mAb with DAB staining.

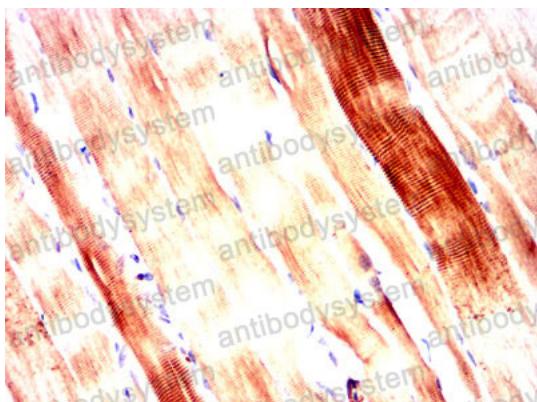

Immunohistochemical

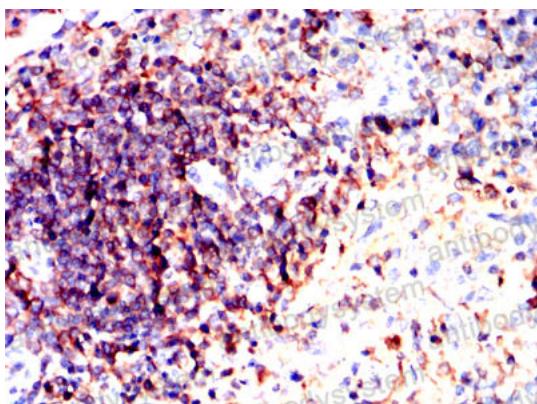

Immunohistochemical

Immunohistochemical analysis of paraffin-embedded Rat myocardium tissues using ATM mouse mAb with DAB staining.

Immunohistochemical analysis of paraffin-embedded Rabbit spleen tissues using ATM mouse mAb with DAB staining.